Specifications

CS-FNX series

Model	CS150FNX	CS120FNX				
Maximum speed (rpm)	150,000	120,000				
Maximum RCF (xg)	1,050,000 (S140AT rotor)	771,000 (S140AT rotor)				
Speed control accuracy (rpm)	+/-50 (5,000 to max. speed)					
Max. capacity (nominal)	30ml x 6 tubes (S50A rotor)					
Max. capacity (nominal) (Swinging bucket rotor)	7ml x 4 tubes (S50ST rotor)					
Timer	1 min. to 99 hrs. and 59 min. with HOLD and RTC (real-time control) function					
Vacuum system	Oil rotary vacuum pump and oil diffusion pump					
Rotor temp. control range (°C)	0 to +40 (1 degree increment)					
Rotor Cooling method	Thermo-module cooling system (HFC free)					
Screen display	Colour LCD (touch-sensitive)					
Rotor setting method	Self-locking rotor system (no special tool required)					
Noise level	45 dB(A) measured 1 meter from front					
Max. heat radiation into room (kW)	0.7 or less					
Drive unit warranty	5 years					
Dimensions (mm)	440 (W) x 520 (D) x 910 (H), height to table 790					
Floor area (m²)	0.53 (750 x 700 mm)					
Weight (kg)	105					
Power requirements		5, 15A (50/60Hz), single phase /- 10%, 8A (50/60Hz), single phase				
Standard	IEC61010-1 and IEC	61010-2-020, EN61326, CE marking qualified				

CS150NX

Model	CS150NX
Maximum speed (rpm)	150,000
Maximum RCF (xg)	1,050,000 (S140AT rotor)
Speed control accuracy (rpm)	+/-50 (5,000 to max. speed)
Max. capacity (nominal)	30ml x 6 tubes (S50A rotor)
Max. capacity (nominal) (Swinging bucket rotor)	7ml x 4 tubes (S50ST rotor)
Timer	1 min. to 99 hrs. and 59 min. with HOLD and RTC (real-time control) function
Vacuum system	Oil rotary vacuum pump and oil diffusion pump
Rotor temp. control range (°C)	0 to +40 (1 degree increment)
Rotor Cooling method	Thermo-module cooling system (HFC free)
Screen display	Colour LCD (touch-sensitive)
Rotor setting method	Self-locking rotor system (no special tool required)
Noise level	45 dB(A) measured 1 meter from front
Max. heat radiation into room (kW)	0.7 or less
Drive unit warranty	5 years
Dimensions (mm)	590 (W) x 582 (D) x 408 (H)
Floor area (m ²)	0.50 (850 x 590 mm)
Weight (kg)	97
Power requirements	AC 110 or 120V ±10%, 15A (50/60Hz), single phase AC 208, 220, 230, 240 V +/- 10%, 8A (50/60Hz), single phase
Standard	IEC61010-1 and IEC61010-2-020, EN61326, CE marking qualified



Koki Holdings Co., Ltd. is certicated with ISO 14001, Environmental management systems by former name, Hitachi Koki Co., Ltd.

Life-Science Instruments Division of Koki Holdings Co., Ltd. is certificated with ISO 9001, Quality management system by former name, Hitachi Koki Co., Ltd.

For safety and proper use of your machine, carefully read and follow the instruction manual.

ISO 9001

Certification

BUREAU VERITAS

All specifications are subject to change without advance notice. For further information, please contact us.

Ŵ

UŘAS

IANAGEMENT SYSTEMS

008



CAUTION:

Koki Holdings Co., Ltd. Life-Science Instruments Division Shinagawa Intercity Tower A, 2-15-1 Konan, Minato-ku, Tokyo 108-6018 Japan Tel : +81-3-6738-0863 Fax : +81-3-6738-0861 https://www.himac-science.com

We change our brand, himac from Hitachi, effective on June 2018. And change our company name. "Koki Holdings Co., Ltd." from "Hitachi Koki Co., Ltd." on June 2018.

HLI-2E 2018.06 Printed in Japan



CS150NX

CS-FNX SERIES

Micro Ultracentrifuges





"Floor-standing model" or "Tabletop model" ?





CS150NX

Micro Ultracentrifuges Floor-standing

CS150FNX CS120FNX

himac ultracentrifuges with one of top performance in speed and RCF. Solutions applied in life science to separate and purify protein, lipoprotein, cell organelle, uncleic acid, DNA, and RNA. And also carbon nano tubes and other nano-size particles in material science. For pretreatment of proteomic analysis.

- Maximum speed 150,000 rpm for CS150FNX and CS150NX with optional rotor S150AT
- Maximum RCF 1,050,000 xg for CS150FNX and CS150NX with optional rotor S140AT
- Quiet low noise 45 dB(A)
- Powerful imbalance protection
- Quick startup
- Program and step-mode operation







Biosafety Option

In order to prevent from exhausting bio-hazardous sample to a room, a micro filter can be assembled in a vacuum line at option.

Tabletop **CS150NX**

- Administration system with user lockout
- Self-locking rotor system
- Data communication
- Optional rotor with large capacity
- Optional software

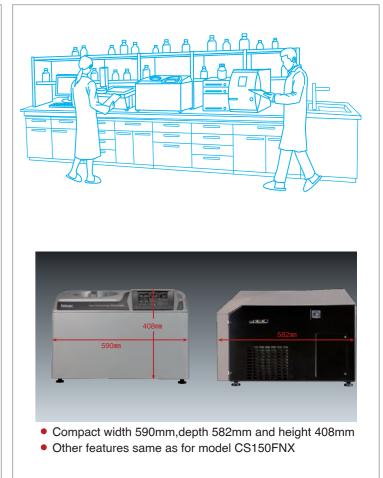

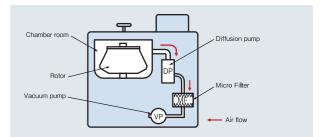

Features & Functions

One of the world's top performance: Speed 150,000 rpm and RCF 1,050,000 xg by CS150FNX and CS150NX*


himac state-of-art technologies has realized the greatest performance above and quiet operating sound low noise 45dB(A) in compact body*. Now you can separate and purify protein, lipoprotein, cell organelle, DNA and RNA, also carbon nano tubes and other nano-sized particles, more efficiently than ever.

The Most Reliable Drive Unit

CS150FNX and CS150NX take only 90 seconds to reach the world's fastest speed 150,000rpm by himac original vacuumed drive unit. Because of its high reliability, we offer 5 years warranty to the drive unit.

Powerful Imbalance Protection

Samples need to be balanced within 5mm by visual check only. Non-contact imbalance sensor always monitors vibration of the rotor and drive shaft. In case of unusual vibration, the sensor activates and stops operation immediately.

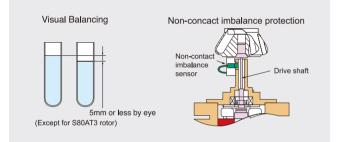

Quick Startup

It takes just 8 seconds (minimum) to be ready for operation since power switch is turned on. Initial screen is displayed on the LCD panel during booting up the system. Reaching time to high-vacuum status is reduced to half by optimized control of the oil diffusion pump. (comparing with our previous model.)



Environment-friendly Design

CS-FNX series and CS150NX are designed according to our "Design for Environment assessment system" and designated as Eco-Products.

60% stand-by electricity and maximum 10% operating power consumption are reduced, compared with our former models. Also over 80% materials of the product are reusable or recyclable.


Easy Operation with Touch-Sensitive LCD

Color touch-sensitive LCD and GUI (graphic user interface) with high contrast against back screen in black color enable users easily to operate the system or select various menus and functions by touching the icon on the display. Of course, operation status is identified at a glance by intelligible screen design.



Quiet Low Noise 45 dB(A)

Quiet operating sound is realized by new rigid drive system and sound absorbing structure designed with CAE analysis. It does not disturb you in laboratory or office.

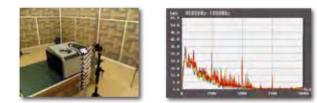

Easy Timer Setting and Actual Run Timer

himac original RTC (real-time control) function makes timer setting easy (PAT.). What you need to do is simply set start time or finish time with running time. It is easier than conventional delay time setting timer.

himac original Actual Run Timer starts when the set speed is attained (PAT.) and excludes acceleration time from the set time. It helps to precisely control the net run-time and obtain high-reproducibility separation. Of course, conventional run timer is selectable.

Administration system with Use Lockout

You can register users to the system with different access levels and password. The system can be protected by password to enter the system. This function is very useful when you need to control users who can access the system.

Data Communication

USB port is equipped as standard specification. (CS150FNX and CS150NX only) The system records up to 100 operating histories in its memory (PAT.). So operating data can be output in CSV format through the USB port.

Ethernet is a registered trademark of Xerox Corporation in the United States and other countries.

himac LogManager for Windows® (Network Edition) - optional log management software

himac LogManager for Windows® (Network Edition) is useful and convenient log management software of operating histories of CS150FNX and CS150NX. Maximum 16 units can be registered to the software. As data communication between the centrifuge and PC is done via Local Area Network. You can establish the network configuration and relocation easily with commercially available LAN devices. Optional LAN board is required to CS150FNX and CS150NX. It means you can manage the operating log at the different location from the installation site of the centrifuges.

This software is real-time log management software, recording period is selectable from 10 seconds to 5 minutes interval. Of course, the software supports USA FDA 21 CFR Part 11, following functions are available;

- Digital Signature
- Audit Trail
- Encrypted Data Files
- Optional for CS150FNX and CS150NX, and not for CS120NX.
- Required Operating System : Windows 7[®] Professional or Windows 8[®] Professional

Optional LAN board is also available. (CS150FNX and CS150NX only) You can link CS150FNX and/or CS150NX to your PC through LAN by Ethernet. Optional software "himac LogManager (network edition)" can be installed in your PC and operating history and data of CS150FNX and/or CS150NX can be easily managed by your PC.

COLUMN 1

Rotors

Self-Locking Rotor System

himac original "self-locking rotor system" is simple system to set the rotor onto the drive shaft. You just simply place the rotor on the drive shaft and that's all. The rotor is locked by centrifugal force automatically. You do not need either to fix the rotor by screw or push the button to lock or unlock the rotor like other models available in the markets.



Large nominal capacity 7 ml x 4 tubes

for ultracentrifuge by 7ml tubes.

Top loading System

Easy to secure a swing bucket to set

onto rotor by top loading to avoid its

50,000

253,000

77

7.0 x 4 7.0 x 4

9124401K

wrong setting.

Larger Capacity Rotors

S50A



Maximum nominal capacity is 30ml x 6 tubes, it is one of the largest capacity rotor for micro ultracentrifuge. Also it is easy to open and close the cover by T-shape handle. Newly introduced 30ml tubes series and himac original S-cap can be used for S50A. **S50ST**

Max. speed (rpr Max. RCF (xg)

Capacity by standard tube (ml x tubes)*1

Max. capacity by optional tube (ml x tubes)

K-facto

Part No

Max. speed (rpm)	50,000
Max. RCF (xg)	210,000
K-factor	61
Capacity by standard tube (ml x tubes)*1	25.0 x 6
Max. capacity by optional tube (ml x tubes)	30.0 x 6
Part No.	9124391K

Note : *1. Nominal capacity by standard tube included with the rotor

S-cap Series

himac original S-cap series offer simple and easy operation of the tube cap. This unique system requires just three components and three steps to assemble the cap to the tube.

- No consumables, such as O-ring
- Easy cleaning

Three assmebling steps

- 1. Fill sample into the tube up to 80% volume and insert AL Stem into the tube by the tool.
- 2. Assemble Ti Ring on the tube and tighten AL Stem and Ti Ring by the tool.
- 3. Fill up the sample into the tube by an injector and close stem by Setscrew.

Exclusive 1.5ml microtube (conical bottom) x 12pcs 201,000 x g has been realized !



S55A2 Protein and nucleic acid etc. are quickly separated and the pretreatment of proteomic analysis is accelerated.

You can use this exclusive microtube up to 201,000 x g without any adaptor in S55A2

1.5ml microtube (applicable to high x g)

rotor.



(300pcs/box) For further details, please refer to himac application No. 127 available on our web site. (http://www.himac-science.com)

(P/N S308420B)

Rotor Type	Model	Part No.	Max. speed	Max. RCF (×g)	K-factor	Capacity by standard tube (ml x tubes)*1	Max Capacity by optional tube (ml x tubes)
Fixed angle	S150AT	9124211K	150,000	901,000	6	1.0×8	*2 2.0×8
	S140AT	9124351K	140,000	1,050,000	5	1.0×10	*2 2.0×10
	S120AT2	9124181K	120,000	650,000	8	1.0×10	*2 2.0×10
	S120AT3	9124191K	120,000	650,000	8	0.5×14	0.5×14
	S110AT	9124411K	110,000	691,000	15	4.0×8	*2 5.0×8
	S100AT3		100,000	436,000	7	0.23×20	0.23×20
	S100AT4	9124241K	100,000	541,000	16	3.0×6	*2 4.0×6
	S80AT2	040405414		358,000	14	0.5×20	0.5×30
	300A12	9124251K	80,000	304,000	17	0.5×10	0.5 ~ 30
	S80AT3	9124331K	80,000	415,000	23	6.0×8	*2 8.0×8
	S70AT	9124261K	70,000	307,000	31	0.5×20	0.5×20
	S58A	9124381K	58,000	289,000	50	10.0×8	*2 13.5×8
	S55A2	9124371K	55,000	201,000	40	1.5×12	*4 1.5×12
	S50A	9124391K	50,000	210,000	61	25.0×6	*1 30.0×6
	S140NT	9124223K	140,000	752,000	6	2.0×8	*3 2.0×8
Neo angle	S120NT	9124303K	120,000	586,000	9	2.0×8	*3 2.0×8
	S100NT	9124313K	100,000	479,000	12	4.0×8	*3 4.0×8
	S55S	9124281K	55,000	259,000	44	2.2×4	2.2×4
Swinging bucket	S52ST	9124341K	52,000	276,000	79	5.0×4	5.0×4
	S50ST	9124401K	50,000	253,000	77	7.0×4	*1 7.0×4
Vertical	S120VT	9124291K	120,000	501,000	8	2.0×8	*3 2.0×8

Notes 1. Nominal capacity by standard tube included with the rotor. 2. Capacity indicated is when using a seal tube. Please place a separate order for STF2 tube sealer, tube rack and space cap. 3. Space cap is standard accessory.

4. Third-party's 1.5ml microtubes cannot be used. 5. ● indicates a made-to-order rotor. 6. ▲ indicates sales discontinued.

Tubes for Rotor

		Tube licates standard tube for rotor	Tube Nominal	Tube s	specificat		Tube supply		ccessories neede e tubes (order se	
Rotor	Part No.	Name	capacity	Size (mm)	Material	Max. speed	unit	STF3	Adapter	Cap
			(ml)	, <i>,</i> ,		(rpm)	(pcs/box)	Tube rack	Auapter	Cap
S150AT	S300535A	1PC tube	1.0	φ 11 x 34	PC	150,000	100			
	S300539A	1.5PA seal tube*	1.5	φ 11 x 25	PA	150,000	50	Tube rack (G2)		Space cap (A4) x 8
	S300540A	2PA seal tube*	2.0	φ 11 x 32	PA	150,000	50	Tube rack (G)***		Space cap (A4) x 8
0440.47	S309652A	1PA Tube	1.0	φ 11 x 34	PA	86,000	100			
S140AT	S300535A	1PC tube	1.0	φ 11 x 34	PC	140,000	100	Tube male (OO)		0
	S300539A	1.5PA seal tube*	1.5	φ 11 x 25	PA	140,000	50	Tube rack (G2)		Space cap (A4) x 10
	S300540A	2PA seal tube*	2.0	φ 11 x 32	PA	140,000	50	Tube rack (G) ***		Space cap (A4) x 10
0100 470	S309652A S300535A	1PA Tube	1.0	φ 11 x 34	PA	86,000	100			
S120AT2		1PC tube • 1.5PA seal tube*	1.0	ϕ 11 x 34	PC PA	120,000	100 50	Tube rack (G2)		Space cap (A4) x 10
	S300539A S300540A	2PA seal tube*	1.5 2.0	φ 11 x 25 φ 11 x 32	PA	120,000	50	Tube rack (G2)		Space cap (A4) x 10 Space cap (A4) x 10
	S309652A	1PA Tube	1.0	ϕ 11 x 32 ϕ 11 x 34	PA	86,000	100	Tube fack (G)		
S120AT3	S309652A S300533A	0.5 PC tube	0.5	ϕ 11 x 34 ϕ 8 x 34	PC	120,000	100			
0120410	S300534A	0.5 PA tube	0.5	φ 8 x 34	PA	80,000	100			
S110AT	S404332A	4PC tube	4.0	ϕ 13 x 6	PC	110,000	50			
SHUAI	S303351A	4PA thick-walled tube	4.0	φισχο	PA	80,000	50			
	345319A	5PA seal tube	5.0	φ 13 x 53	PA	110,000	50	Tube rack (B) ***		B4 space cap x 8
	S308892A	1.5ml micro tube**		ϕ 10.9 x 40	PP	48,000	300	Tube Tack (D)	1.5A5 adapter x 4	D4 Space cap x 0
S100AT3	339133A	0.23 PC tube	0.23	ϕ 10.9 x 40 ϕ 7 x 21	PC	100,000	100			
STODAT3	339133A 339675A	0.23 PC tube	0.23	ϕ 7 x 21 ϕ 7 x 21	PC	80,000	100			
S100AT4	S308892A	1.5 micro tube**	1.5	ϕ 11 x 42	PA PP	50,000	300		1.5A3 adapter x 6	
5100A14	S308892A S301599A	3PC tube	3.0	$\frac{\phi}{\phi} \frac{11 \times 42}{13 \times 51}$	PP PC	,	50		1.575 auapter x 6	
				,		100,000		Tubo rock (DO) ***		B4 cross car v 2
	S303694A	4PA seal tube*	4.0	ϕ 13 x 38	PA	100,000	50	Tube rack (B3) ***		B4 space cap x 6
0400.0=4	S303350A	3PA thick-walled tube	3.0	ϕ 13 x 51	PA	85,000	50			
S100AT6	S308892A	1.5 micro tube**	1.5	φ 11 x 42	PP	47,000	300		1.5A5 adapter x 4	
Sales	S404332A	4PC tube	4.0	φ 13 x 60	PC	100,000	50			
discontinued)	L	5PA seal tube*	5.0	φ 13 x 53	PA	100,000	50	Tube rack (B) ***		B4 space cap x 8
	S303351A	4PA thick-walled tube	4.0	φ 13 x 60	PA	80,000	50			
S80AT2	S300533A	0.5 PC tube	0.5	φ 8 x 34	PC	80,000	100			
	S300534A	0.5 PA tube	0.5	φ 8 x 34	PA	80,000	100			
S80AT3	S308892A	1.5 micro tube**	1.5	φ 11 x 42	PP	45,000	300		1.5A8 adapter x 8	
	S408829A	6PC thick-walled tube	6.0	ϕ 16 x 59	PC	80,000	50			
	S408833A	6PA thick-walled tube	6.0	ϕ 16 x 59	PA	60,000	50			
	S408831A	8PA seal tube*	8.0	φ 16 x 53	PA	80,000	50	Tube rack (C2)		Space cap (C5) x 8
S70AT	S302371A	0.5 PC tube (B)	0.5	φ7 x 56	PC	70,000	100			
S58A	345320A	12PA seal tube*		φ 16 x 78	PA		50	Tube rack(C)***		C2 space cap x 8
	329606A	12PA tube		φ 16 x 76	PA	E9 000	50			C-TI cap ass'y x 8
	325470A	12PE tube		φ16 x 76	PE	58,000	50			C-TI cap ass'y x 8
	325751A	12PC tube	12.0	φ 16 x 76	PC		50			C-TI cap ass'y x 8
	0000404	1000 to be		(1070	00	30,000				C-TI cap ass'y x 8
	326010A	12SS tube		φ 16 x 76	SS	35,000	2			
	338455A	10PC thick-walled tube		φ 16 x 76	PC	58,000	10			
	338456A	10PA thick-walled tube		φ 16 x 76	PA	50,000	10			
	S303922A	10PC bottle (B) ass'y.	10.0	φ 16 x 82	PC	58,000	6			
	325952A	10PC bottle ass'y.		φ 16 x 82	PC	55,000	10			
	334105A	10PA bottle ass'y.		φ 16 x 82	PA	50,000	10			
	329445A	6.5PA tube		φ 13 x 64	PA	38,000	50		6.5A12 adapter	B-AL cap ass'y x 8
	326415A	6.5PC tube	6.5	φ 13 x 64	PC	38,000	50		6.5A12 adapter	B-AL cap ass'y x 8
	332245A	5PA tube	5.0	φ 13 x 52	PA	38,000	50		5A12 adapter	B-AL cap ass'y x 8
	S303141A	4.7PC thick-walled tube	4.7	ϕ 13 x 64	PC	38,000	50		6.5A12 adapter	
	329561A	4PA tube	4.0	ϕ 10 x 64 ϕ 11 x 60	PA	38,000	50		4A12 adapter	A3-AL cap ass'y
	336708A	3PA tube		ϕ 11 x 47	PA	38,000	50		3A12 adapter	A3-AL cap ass'y
	S303271A		3.0	ϕ 11 x 47	PET	38,000	50		3A12 adapter	A3-AL cap ass'y
	S301599A		0.0	ϕ 13 x 51	PC	38,000	50		5A12 adapter	10 11 0up ass y
	S304296A		0.9	ϕ 13 x 51 ϕ 8 x 50	PC	42,000	50		2A12 adapter	
S55A2	S304290A S308892A	1.5 micro tube	1.5	ϕ 11 x 42	PP	55,000	300		Linz adapter	
555A2 S50A	S308892A S309154A	30PA tube ass'y	30.0	$\frac{\phi}{\phi} \frac{11 \times 42}{26 \times 72}$	PP PA	50,000	50			E3-AL cap ass'y x 6
	00004404	05DO thigk well at the state	05.0			50.000	10			S-40AL cap ass'y x
		25PC thick-walled tube ass'y●	25.0	φ 26 x 76	PC	50,000	10			
		25PA thick-walled tube ass'y	0.5.4	1.65.53	PA	30,000	10			
		20PC bottle(C) ass'y	20.0	φ 26 x 80	PC	50,000	6	-		
		2PA seal tube*	2.0	φ 11 x 32	PA	140,000	50	Tube rack (G) ***		Comes with rotor
	S300540A			φ 11 x 32	PA	120,000	50	Tube rack (G) ***		Comes with rotor
S120NT	S300540A	2PA seal tube*	2.0	1.1.7	PA	100,000	50	Tube rack (B3)***		Comes with rotor
S140NT S120NT S100NT	S300540A S303694A	2PA seal tube* ● 4PA seal tube* ●	4.0	φ 13 x 38				Lubo rook (C) ***		
S120NT S100NT S120VT	S300540A S303694A S300540A	2PA seal tube* 4PA seal tube* 2PA seal tube*	4.0 2.0	φ 11 x 32	PA	120,000	50	Tube rack (G)***		Comes with rotor
S120NT S100NT S120VT	S300540A S303694A S300540A S300535A	2PA seal tube* 4PA seal tube* 2PA seal tube* 1PC tube	4.0 2.0 1.0	φ 11 x 32 φ 11 x 34	PA PC	55,000	100			Comes with rotor
S120NT S100NT S120VT	S300540A S303694A S300540A	2PA seal tube* 4PA seal tube* 2PA seal tube*	4.0 2.0	φ 11 x 32	PA					Comes with rotor
S120NT S100NT S120VT	S300540A S303694A S300540A S300535A	2PA seal tube* 4PA seal tube* 2PA seal tube* 1PC tube	4.0 2.0 1.0	φ 11 x 32 φ 11 x 34	PA PC	55,000	100			Comes with rotor
S120NT S100NT S120VT S55S	S300540A S303694A S300540A S300535A S300536A	2PA seal tube* 4PA seal tube* 2PA seal tube* 1PC tube 2.2PA tube	4.0 2.0 1.0 2.2	φ 11 x 32 φ 11 x 34 φ 11 x 34	PA PC PA	55,000 55,000	100 50			Comes with rotor
S120NT S100NT S120VT S55S	S300540A S303694A S300540A S300535A S300536A S300538A 332245A	2PA seal tube* 4PA seal tube* 2PA seal tube* 1PC tube 2.2PA tube 2.2PET tube	4.0 2.0 1.0 2.2 2.2	φ 11 x 32 φ 11 x 34 φ 11 x 34 φ 11 x 34	PA PC PA PET	55,000 55,000 55,000	100 50 50			Comes with rotor
S120NT	S300540A S303694A S300540A S300535A S300536A S300538A 332245A	2PA seal tube* 4PA seal tube* 2PA seal tube* 1PC tube 2.2PA tube 2.2PET tube 5PA tube	4.0 2.0 1.0 2.2 2.2 4.9	$ \phi 11 x 32 \phi 11 x 34 \phi 11 x 34 \phi 11 x 34 \phi 13 x 52 $	PA PC PA PET PA	55,000 55,000 55,000 52,000	100 50 50 50			Comes with rotor
S120NT S100NT S120VT S55S	S300540A S303694A S300540A S300535A S300536A S300538A 332245A S303273A	2PA seal tube* 4PA seal tube* 2PA seal tube* 1PC tube 2.2PA tube 2.2PET tube 5PA tube 5PET tube 3PC tube	4.0 2.0 1.0 2.2 2.2 4.9 5.3	$ \phi 11 \times 32 \\ \phi 11 \times 34 \\ \phi 11 \times 34 \\ \phi 11 \times 34 \\ \phi 13 \times 52 \\ \phi 13 \times 52 \\ \phi 13 \times 52 \\ \end{cases} $	PA PC PA PET PA PET	55,000 55,000 55,000 52,000 52,000	100 50 50 50 50		2S5 adapter x 4	Comes with rotor
S120NT S100NT S120VT S55S	S300540A S303694A S300540A S300535A S300536A S300538A 332245A S303273A S301599A S304296A	2PA seal tube* 4PA seal tube* 2PA seal tube* 1PC tube 2.2PA tube 2.2PET tube 5PA tube 5PET tube 3PC tube 0.9PC thick-walled tube	4.0 2.0 1.0 2.2 2.2 4.9 5.3 3.4	$ \phi 11 \times 32 \phi 11 \times 34 \phi 11 \times 34 \phi 11 \times 34 \phi 13 \times 52 \phi 13 \times 52 \phi 13 \times 51 \phi 8 \times 50 $	PA PC PA PET PA PET PC PC	55,000 55,000 55,000 52,000 52,000 52,000 43,000	100 50 50 50 50 50 50 50		2S5 adapter x 4	Comes with rotor
S120NT S100NT S120VT S55S S52ST	S300540A S303694A S300540A S300535A S300536A S300538A 332245A S303273A S301599A	2PA seal tube* 4PA seal tube* 2PA seal tube* PC tube 2.2PA tube 2.2PET tube 5PA tube 5PET tube 3PC tube 0.9PC thick-walled tube 7PA tube(B)	4.0 2.0 1.0 2.2 2.2 4.9 5.3 3.4 0.9	$ \phi 11 \times 32 \phi 11 \times 34 \phi 11 \times 34 \phi 11 \times 34 \phi 13 \times 52 \phi 13 \times 52 \phi 13 \times 51 $	PA PC PA PET PA PET PC	55,000 55,000 55,000 52,000 52,000 52,000	100 50 50 50 50 50 50		2S5 adapter x 4	Comes with rotor

Notes: 1. When using tubes marked with *, optional accessories are required (STF3 tube sealer, tube rack and space cap). 2. Tubes marked with ** are commercially available too. When using these items, 1.5A3, 1.5A5 or 1.5A8 adapters are necessary; therefore, please order separately. 3. ***Also applicable with Multi rack, Part No.S201778K. 4.

indicates tube included with its rotor as its standard accessory.

Three Components S-cap tool